


Perineal & Episiotomy Repair Trainer

Product No: 60450


A three stage training system for teaching episiotomy, perineal suturing techniques, and repair of episiotomy and second degree tears.

Developed in collaboration with:


- Professor Christine Kettle, Professor of Women's Health, University Hospital of North Staffordshire and Staffordshire University (UK).
- Professor Khaled Ismail, Senior Lecturer and Consultant Obstetrician and Gynaecologist, Keele Medical School and University Hospital of North Staffordshire (UK).

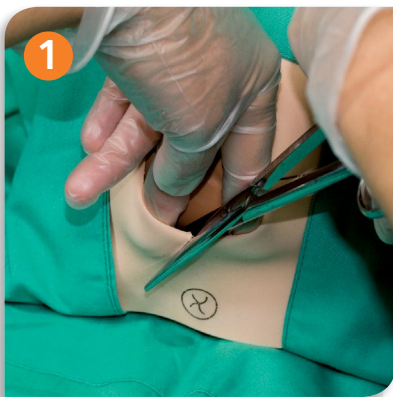

1 Episiotomy

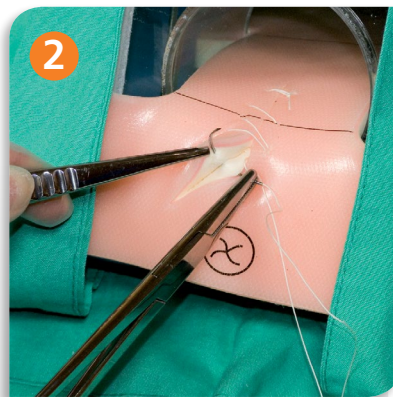

2 Perineal Repair Techniques

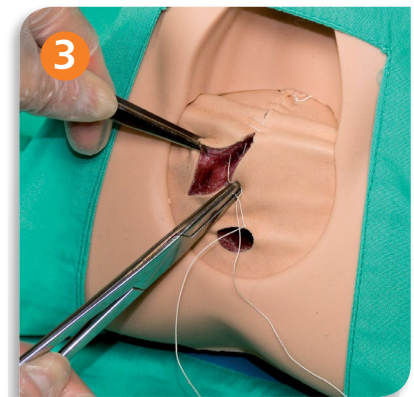

3 Perineal Repair Procedures

Learning has never been more lifelike

limbsandthings.com

Stage 1 Episiotomy

A simplified model for teaching safe methods of episiotomy.

- Identification of fontanelles
- Handling of a stretched perineum
- Reducing the risk of damaging the baby or mother
- Infiltration of perineum prior to episiotomy
- Performance of medio-lateral (or midline) episiotomy
- Use of instruments



– Realistic representation of tissue to represent a stretched perineum



– Baby head provides a realistic look and feel providing landmarks and appropriate tactile response when inserting fingers to guard the baby's head

Stage 2 Perineal Repair Techniques

The Episiotomy and Perineal Repair Pad prompts trainees to learn to suture on 2 planes, one of which has the spatial challenges which exist when suturing within the vagina. Comprises epidermis, dermal and subdermal layer.

- Suturing in 2 planes: inside vagina and on the perineum
- Practice of suture techniques: continuous, subcuticular, knot tying



– Soft skin with a similar drag and strength to human skin
– Epidermis and dermis have a realistic retention of suture



– Representation of the anus printed on pad for orientation
– Representation of the hymenal remnant printed on pad to assist in assessment of suturing techniques

Stage 3 Perineal Repair Procedures

For the practice of episiotomy and second degree tear repair.

- Tissue layer recognition and handling
- Vaginal mucosal suturing
- Deep muscle suturing
- Subcuticular suturing
- Identification and management of perineal tears
- Enables digital rectal examination before and after repair



– Perineal skin
– Superficial muscle structure and relevant layer for suturing



– Cost effective: all layers can be sutured repeatedly
– Each Procedure Block presents 3 repair sites (left lateral is pre-incised, with both midline and right lateral uncut)

You'll Receive

- **60451** 1 Episiotomy Incision Pad (x6)
- **60452** 1 Episiotomy & Perineal Repair Pad (x2)
- **60453** 1 Episiotomy & Perineal Repair Block
- **60454** 1 Perineal Repair Perineum
- **60455** 1 Perineal Repair Techniques Jig
- **60456** 1 Perineal Repair Base & Clip
- **60457** 1 Keele & Staffs Baby Head
- 2 Bench top clamps
- 1 Carry case

Packed weight (kg) & dimensions (mm)

4.4 kg

460 l

310 w

180 h

For demonstration or to order



Limbs & Things Ltd. Sussex Street, St Phillips, Bristol, BS2 0RA, UK

T: 00 44 (0)117 311 0500 F: 00 44 (0)117 311 0501 E: sales@limbsandthings.com

limbsandthings.com

V.2 Mar 2018